

# PROFESSIONAL REHAB- AND PHYSIOTHERAPY-EQUIPMENT

## The Vacuumcushion program for clinic, doctors office and rehabilitation

4FORM® VACUUMCUSHION – THE GENTLE PATIENT BEDDING


<p><b>4FORM®</b> vacuum cushion 30 x 40 cm # 461 3040 e.g. for positioning of hand, head, foot or forearm</p>	<p><b>4FORM®</b> vacuum cushion 30 x 70 cm # 462 3070 e.g. for positioning arm and leg</p>	<p><b>4FORM®</b> vacuum cushion 40 x 70 cm # 463 4070 e.g. for positioning arm and leg or in the incubator</p>	<p><b>4FORM®</b> vacuum cushion 60 x 80 cm # 464 6080 e.g. for the back position- ing or for infants</p>	<p><b>4FORM®</b> vacuum cushion 75 x 100 cm # 465 75100 e.g. for the back positioning or for toddlers</p>	<p><b>4FORM®</b> vacuum cushion 75 x 200 cm # 466 75200 for positioning of the complete body</p>
---	--	--	--	---	--

SPRECHEN SIE UNS GERNE AUF SONDERANFERTIGUNGEN AN!

The individual and ergonomic body adaptation of **4FORM** vacuum cushions serves physically handicapped persons a safe support in different body positions. Their individual formability is achieved by Styrofoam balls, which stabilize by vacuum to adapted constellations. Useable for X-ray, CT or MRI treatment. All vacuum cushions are fitted out with an easy managing two-part-CPC valve.

They can also be use as styrofoam-cushion when the valve is open, and when you blow the cushion up, you can use it as air cushion.

All **4FORM®** vacuum cushions are made to be non-slip and dimensionally stable from a special soft Polyurethan-foil. They are washable and disinfectable, thats why they are hygienic.

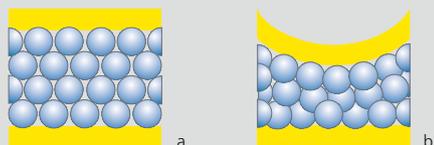

### Styrofoam balls constellation

- without vacuum
- with vacuum



### ACCESSORY FOR 4FORM VACUUMCUSHIONS



Repairkit # 3490124  
For repairing small tears and holes in all 4FORM® vacuumcushions (PU-foil, blue)



Plastic vacuum pump  
# 5370 length 45 cm, red  
# 5355 anti-magnetic brass rod,  
length 45 cm, red



Plastic products for health and emergency services

